



# SIDEKICK SHARPENER

*Sharpening innovation at your fingertips*



## FACT

Scalers and curettes can dull over time as a result of routine, every day use. Keeping instruments sharp can improve deposit removal and decrease the potential for hand fatigue\*.

## CHALLENGE

Achieving a fine cutting edge without compromising the original shape of the instrument blade can be challenging and time-consuming.

## SOLUTION



Hu-Friedy's Sidekick is a cordless, lightweight sharpening device designed to provide convenient and consistent sharpening results and accurate angles every time.

## POINTS OF PERFORMANCE

- Easy-to-use guide channels and vertical backstops allow simple positioning of scaler and curettes to provide consistent sharpening results every time.
- Provides consistent sharpening results aiding in instrument longevity.
- Cordless, lightweight, portable sharpening...anywhere, anytime.
- Reciprocating stone does the work so you don't have to
- Maintains sharp blade and accurate angles




# SIDEKICK SHARPENER

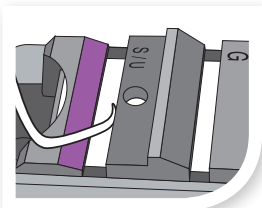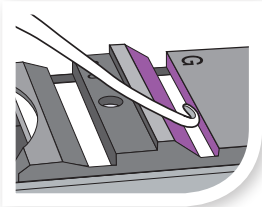


## SHARPENING MADE EASY!

- Easy-to-read letters indicate Graceys (**G**) or Sickles and Universals (**S/U**)
- Position the instrument in the specific channel, with the terminal shank resting on the incline of the channel
- Position the back of the instrument along the backstop of the guideplate
- Turn the unit on, and glide the instrument within the channel from side to side
- Repeat 2 - 3 times or until blade is sharp


## ACHIEVE ACCURATE ANGLES EVERY TIME

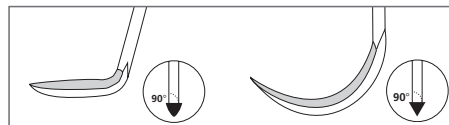

SICKLE SCALARS

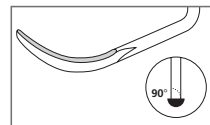

UNIVERSAL CURETTES

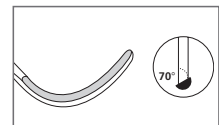

GRACEY CURETTES

## ADDITIONAL ITEMS

### REPLACEMENT PARTS

Stone Clip	SDKCLIP
Guideplate Screw	SDKGPSCREW
Guideplate	SDKPLATE
Ceramic Stone	SDKSTONEC
Sharpen-EZ™ Stone Oil	SSO
Plastic Test Sticks	PTS
Magnifying Lens	LENS
Large Parts Box	IMS-1273

### SIDEKICK SHARPENING KIT: SDKKIT

Kit includes: Sidekick Unit, User Manual, Ceramic Sharpening Stone, Plastic Test Stick, Magnifying Lens and 2 AA Batteries.

Visit us online at [Hu-Friedy.eu](http://Hu-Friedy.eu)

Manufacturer: Hu-Friedy Mfg. Co., LLC | 3232 N. Rockwell Street | Chicago, IL 60618 | USA  
Customer Care Department: Free Phone 00800 48 37 43 39  
Free Fax 00800 48 37 43 40 | Email: [info@hu-friedy.eu](mailto:info@hu-friedy.eu)

